




**USER MANUAL – EN**

**IN 24985 Collagen lamp inSPORTline Coladan**

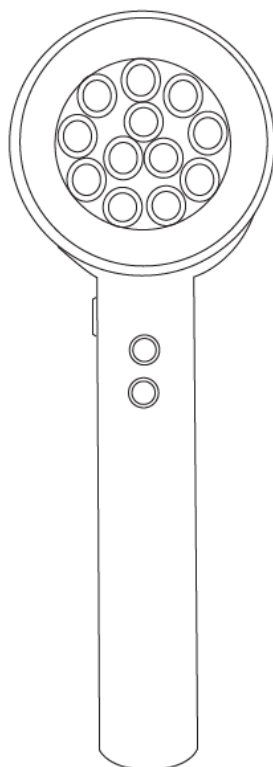

**CONTENTS**

SAFETY INSTRUCTIONS..... 3  
PRODUCT DESCRIPTION ..... 3  
USE ..... 4  
CLEANING AND MAINTENANCE ..... 4  
SPECIFICATIONS..... 4  
ENVIRONMENT PROTECTION ..... 5  
TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS..... 5

## SAFETY INSTRUCTIONS

- Read the manual before use and keep it for future reference.
- First, test the device on your forearm to see if you are allergic or sensitive to light. Then wait about an hour and check the location. If there are no signs of allergies or sensitivities, you can use the device on your face.
- Store the device in a dry and shady place.
- Keep the device away from water.
- Do not use the device while charging.
- Keep out of reach of children and pets.
- Before use, clean the skin and apply a moisturizer. Never use cosmetics containing alcohol.
- Do not use on the following areas: eyes, eyelids, irritated or damaged skin, scars, open wounds, neck, itchy areas on the skin, etc.
- Do not use the device if you are sensitive to light.
- Not suitable for pregnant women or nursing mothers.
- Do not use after plastic surgery.
- Do not use while driving, walking or other activities that require your full attention.

## PRODUCT DESCRIPTION



1. LED lights

2. Intensity adjustment

3. Power switch / mode setting

4. USB-C charging port

5. Device status indicator

**LED lights:** 24 pcs, 630 nm on one side, 660 nm on the other

### Intensity adjustment

**Power / Modes:** Hold down the button to turn the device on / off; after start-up, press the button to switch the set of lights

### USB-C for charging

**Indicator light:** white - normal operation, red (flashing) - low battery, white (flashing) - charging, red - fully charged

### PACKAGE CONTENTS

Apparatus; charging cable; manual; glasses

With regular use, the light has a beneficial effect on the skin (reduction of wrinkles, acne and facial skin imperfections) and on the whole body – it helps relieve pain in muscles, joints and bones, and supports the production of collagen in the body.

## USE

Move the device on your forehead from left to right.



Move from chin to jaw.



Move upwards.



Move around the eyes.



- Recommended time of use: max. 3 min in the same place on the face
- Frequency of use: 5 times a week

## CLEANING AND MAINTENANCE

- Use a dry soft cloth for cleaning.
- The device must not come into contact with water, do not use cleaning agents.

## SPECIFICATIONS

Packaging material	Aluminum alloy
Wavelength of red light	Divided into two halves 630 nm and 660 nm
Number of LED lights	24 pcs
Angle of light	25 °
Power supply	Higher wavelength 15W, lower wavelength 12W
USB-C	5V / 1A
Weight	245 g
Battery capacity	5000 mAh
Dimension	190 x 68 x 35 mm

## ENVIRONMENT PROTECTION

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyards.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

## TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

### General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Strakonická street 1151/2c, Prague 150 00, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

"The Buyer who is the End Customer" or simply the "End Customer" is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

"The Buyer who is not the End Customer" is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

### Warranty Conditions

#### Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

#### Batteries

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

#### The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster

- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

### **Warranty Claim Procedure**

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

#### **CZ**

##### **SEVEN SPORT s.r.o.**

Registered Office: Strakonická 1151/2c, Praha 5, 150 00, ČR  
 Headquarters: Dělnická 957, Vítkov, 749 01  
 Warranty & Service: Čermenská 486, Vítkov 749 01

CRN: 26847264  
 VAT ID: CZ26847264

Phone: +420 556 300 970  
 E-mail: eshop@insportline.cz  
 reklamace@insportline.cz  
 servis@insportline.cz

Web: www.inSPORTline.cz

#### **SK**

##### **inSPORTline s.r.o.**

Headquarters, warranty & service center: Električná 6471, Trenčín 911 01, SK

CRN: 36311723  
 VAT ID: SK2020177082

Phone: +421(0)326 526 701  
 E-mail: objednavky@insportline.sk  
 reklamacie@insportline.sk  
 servis@insportline.sk

Web: www.inSPORTline.sk