



**USER MANUAL – EN**

**IN 24987 LED face panel inSPORTline Coladome 900**



# CONTENTS

- SAFETY INSTRUCTIONS..... 3
- PRODUCT DESCRIPTION ..... 4
- ASSEMBLY ..... 4
- USE ..... 4
  - BUTTONS DESCRIPTION ..... 5
- LIGHT THERAPY ..... 7
- MAINTENANCE ..... 7
- SPECIFICATIONS..... 7
- ENVIRONMENT PROTECTION ..... 7
- TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS..... 8

## SAFETY INSTRUCTIONS

- Read the manual before use and keep it for future reference.
- Before use, check the condition of the device and the LED lights.
- Use the device in a well-ventilated room, avoiding places with high temperature, humidity, and direct sunlight.
- Place the device on a solid, level, and stable surface in a well-ventilated room and out of direct sunlight.
- Protect the device against impacts and falls, move the device slowly and carefully.
- Do not use the device together with other heat-emitting devices.
- Keep out of reach of water.
- Do not use if the power cord or plug is damaged.
- Switch off the device after use or during maintenance and storage.
- If the machine makes abnormal noises or smells, stop using the machine immediately and unplug it.
- Use the device only for the purposes described in the manual.
- Do not repair or modify the device yourself.
- The device is not suitable for people suffering from: heart disease, high blood pressure, skin disease, sensitivity to light, etc.
- If you feel sick, stop using the device immediately.
- Keep out of reach of children and animals.
- Dry your hands thoroughly before use.
- Keep away from flammable substances, liquids, and gases.
- Clean the device with a dry soft cloth. Do not use aggressive cleaning agents.
- Do not place other objects on the device.
- Persons with personality, motor or sensory disorders may only use the device under the supervision of an adult and a responsible person or on the advice of a doctor.
- Do not use unapproved accessories.
- Before first use, we recommend testing the device on the forearm if no allergic reaction or sensitivity to light occurs. Check the tested part after 1 hour, if there are no signs of allergies or sensitivities, you can use the device.
- Always wear safety goggles!
- Recommended operating temperature: 0 - 26 ° C
- Recommended storage temperature: 0 - 35 ° C
- Relative humidity: <80%
- Frequency of use: can be used 20-30 minutes a day for one month or after any cosmetic procedure. As part of acne treatment, we recommend using it 3 times a week, for at least 25 minutes

## PRODUCT DESCRIPTION



- 1. Frame
- 2. Control panel
- 3. Power button
- 4. Vent
- 5. LED lights
- 6. Adapter
- 7. Goggles

With regular use, the lights work as a disinfectant and removes unwanted bacteria from the skin.

## ASSEMBLY




1. Attach both side panels to the main panel with the control panel.
  2. After assembling the device, connect the power cord to the device, and then plug the power cord into an electrical outlet.
  3. An audible signal sounds after a successful power on.
  4. If you have connected the device incorrectly, you will not be able to control it.
- Before use, make sure that you do not place the device near other heat-emitting devices or medical devices.
  - Avoid placing the device in places with high humidity, place the device in a well-ventilated room out of direct sunlight.

## USE

### Infrared radiation

Recommended time of use: 20 - 45 min / one session (adjust the length according to each user)


Infrared light can be used for the whole body.

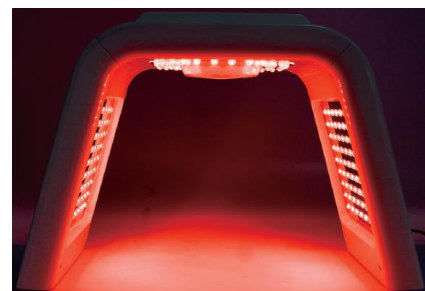


## Ultrasonic cleaning

We recommend using ultrasonic cleansing together with a face mask or moisturizing the face with clean water.

Do not use highly concentrated cosmetic products or essential oils.

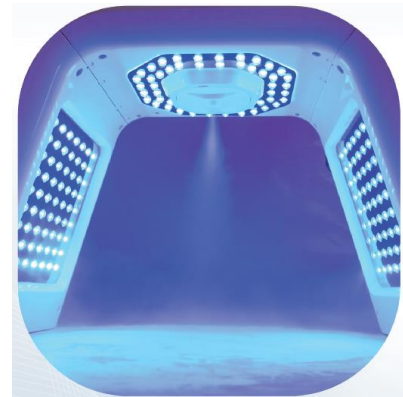

## BUTTONS DESCRIPTION

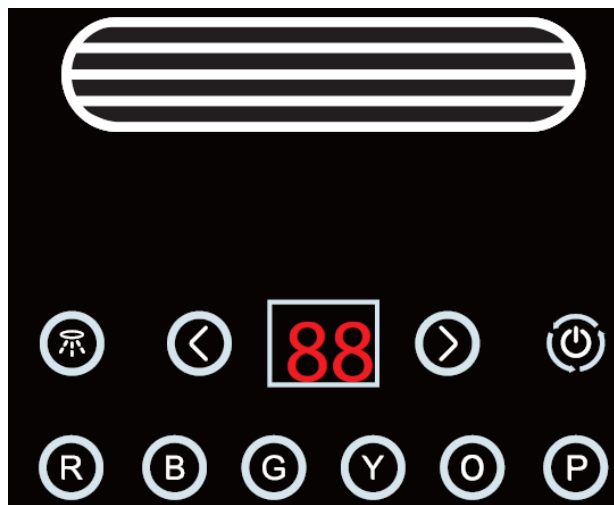

### Power switch / start of automatic function cycling

After connection, the device switches to standby mode. Press the button to start the device and display the operating time of 30 minutes, press it again to switch the device to standby mode.

Hold down the button for 2 - 3 seconds to start the function cycling function.

- Note:**
1. If the cycling function is activated, the colours will change in order (red - orange - yellow - green - turquoise - blue - purple).
  2. The cycling function cannot be started during red, orange, yellow, green, blue and purple light.
  3. If the cycling function is selected, the intensity is set to level 05 and cannot be adjusted.



### Infrared light (wavelength 647 - 760 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Blue light (wavelength 430 - 455 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Green light (wavelength 492 - 550 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Orange light (wavelength 588 - 647 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Yellow light (wavelength 550 - 580 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Purple light (wavelength 360 - 430 nm)



The light intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05).



### Ultrasonic cleaning

The function can be run separately or together with other functions. The intensity can be adjusted in 5 levels 1 (lowest) - 5 (highest). Press the button to adjust the intensity, the set intensity will be shown on the display (01/02/03/04/05). If 00 is displayed, the function is switched off.

### Time +

The default set time is 30 min. After selecting the light function, press the button to increase the time by 5 minutes up to a maximum of 60 minutes. If you exceed the maximum time, the time starts again from 00. The adjustable time range is from 00 - 60 min. If the button is not pressed for 3 seconds, the device starts the set function with the currently set time.

### Time -

The default set time is 30 min. After selecting the light function, press the button to reduce the time by 5 minutes to a minimum of 00 minutes. If you exceed the minimum time, the time starts again from 60. The adjustable time range is from 00 - 60 min. If the button is not pressed for 3 seconds, the device starts the set function with the currently set time.

## LIGHT THERAPY

### Red

Increased blood flow and collagen production. Helps with tissue repair (replacement of old skin cells with new ones) and pore reduction.

### Green

Helps balance the color of the pigment, speed up skin regeneration and reduce oiliness

### Blue

Soothing, anti-inflammatory, destroys skin bacteria

### Orange/Yellow

Energizes the skin, helps reduce wrinkles, balance the color of the pigment and smooth the skin

### Purple

Helps prevent acne

## MAINTENANCE

- Clean the device with a soft, dry cloth.
- Store the device unfolded in a dry and shady place.
- Before any maintenance and handling, the device must first be disconnected from the power supply.
- In the event of any malfunction, contact the complaints / service department.

## SPECIFICATIONS

Made	China
Power source	DC12V, 2A, 50Hz/60Hz
Power supply	Ca. 24W
Dimension	455 x 220 x 315 mm
Weight	1,9 kg ± 10%

## ENVIRONMENT PROTECTION

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyards.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

## **TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS**

### **General Conditions of Warranty and Definition of Terms**

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Strakonická street 1151/2c, Prague 150 00, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

“The Buyer who is the End Customer” or simply the “End Customer” is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

“The Buyer who is not the End Customer” is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

### **Warranty Conditions**

#### **Warranty Period**

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

#### **Batteries**

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

#### **The Warranty does not cover defects resulting from (if applicable):**

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.



## Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

### CZ

#### SEVEN SPORT s.r.o.

Registered Office: Strakonická 1151/2c, Praha 5, 150 00, ČR  
Headquarters: Dělnická 957, Vítkov, 749 01  
Warranty & Service: Čermenská 486, Vítkov 749 01

CRN: 26847264  
VAT ID: CZ26847264

Phone: +420 556 300 970  
E-mail: eshop@insportline.cz  
reklamace@insportline.cz  
servis@insportline.cz

Web: www.inSPORTline.cz

### SK

#### inSPORTline s.r.o.

Headquarters, warranty & service center: Električná 6471, Trenčín 911 01, SK

CRN: 36311723  
VAT ID: SK2020177082

Phone: +421(0)326 526 701  
E-mail: objednavky@insportline.sk  
reklamacie@insportline.sk  
servis@insportline.sk

Web: www.inSPORTline.sk